

STAT PROFILE®
Prime

Exceptional Simplicity in a New Generation Electrolyte Analyser

Introducing an advanced technology electrolyte analyser that has all sensors combined in a single, replaceable, no maintenance card for a simpler, easier analyser.

Complete electrolyte profile
Na, K, Cl, iCa, iMg

**Whole blood, plasma, serum,
urine, and CSF samples**

1-minute stat testing

Optional 10-position autosampler



nova®
biomedical

Easy to Use and Maintain

Nova's advanced technology analyser consists of a no maintenance MicroSensor card that contains all sensors and a ready-to-use, snap-in calibrator cartridge that contains all reagents and a waste container.



Nova MicroSensor cards and calibration cartridges offer a true no maintenance solution for electrolyte testing.

Maximum Readiness and Uptime

Always calibrated and ready

The MicroSensor card is automatically calibrated throughout its 32-day use life and is always ready to deliver an electrolyte profile in just 60 seconds.

Clot Block™ protection

A unique Clot Block sample flow path is designed to protect sensor cartridges from downtime and maintenance that can result from the introduction of a clotted sample.

Safe, maintenance-free waste system

Biohazardous waste is self-contained within the calibrator cartridge, eliminating waste disposal maintenance and potential exposure to biohazardous waste.




Safe, Easy Sampling from Syringes, Tubes, and Cups

A unique safety sample port protects the user from accidental contact with the sample probe and is easily accessed for all sample containers.



Syringes can be docked and then sampled with hands-free operation.

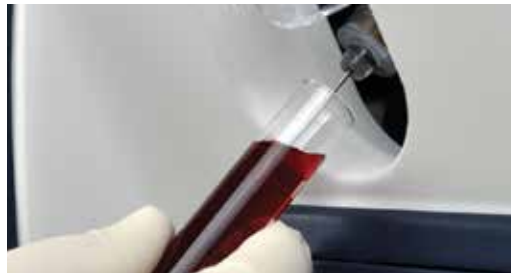


Samples can be aspirated directly from tubes. Sample transfer to a syringe or capillary is eliminated.



A 10-position sample tray accommodates serum, plasma, and urine samples.

Simple, Fast Operation

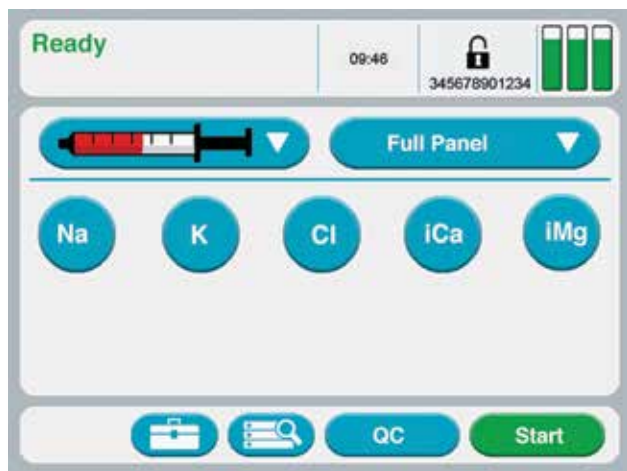

Home Screen

Easy-to-use, high definition, colour touchscreen operation

The 14.2 cm (5.6 in.) colour touchscreen is easily operated through the use of simple and intuitive prompts and requires minimal training.

Three simple steps to operate

1. Press “Start” 2. Scan or enter patient ID 3. Press “Aspirate”

Optional, integrated barcode scanner

An optional 1D/2D barcode scanner, conveniently located within the sample port, eliminates the need for external handheld scanners and allows for fast, error-free entry of operator and patient IDs.

Stat electrolyte results in 60 seconds

Whole blood, serum, or plasma stats are available in just 60 seconds. Whole blood stats can be aspirated directly from syringes or tubes, saving time and sample transfer steps.

New MicroSensor Card Technology



All sensors are contained in one MicroSensor card. Actual size shown.

Accuracy

All Stat Profile Prime tests use proven Nova technology in a miniaturised, sensor-card format. At the end of its 32-day use life, the credit-card sized MicroSensor card is replaced in minutes.

No maintenance

All five electrolyte sensors are built into one small card that is replaced monthly. All individual sensor maintenance procedures, such as membrane replacements, cleaning and polishing procedures, O-ring changes and sensor replacements, are eliminated.

Complete electrolyte profile

Sodium, potassium, ionized calcium, ionized magnesium, and chloride are all required to maintain normal cardiac function. Disruptions in the balance of any of these electrolytes can cause cardiac arrhythmias, reduced cardiac contraction, and complete cardiac arrest. Stat Profile Prime is the only electrolyte system to provide a full assessment of electrolyte balance.



QC statistics and reports are automatically maintained and easily accessed.

True Liquid Quality Control

True liquid QC provides the only reliable test of an analyser

In the United States, federal government regulations (CLIA) are phasing out electronic equivalent quality control (EQC) and are requiring true, liquid-based quality control (QC). Stat Profile Prime combines both true liquid QC and continuous electronic self-monitoring to ensure lab accuracy and uncompromised quality.

QC prompting and QC lockout

Stat Profile Prime provides an on-screen reminder of your customized QC schedule. Out-of-range QC results are automatically flagged and an optional QC lockout feature prevents the reporting of patient results when QC is out of range.

Levey-Jennings charts and QC statistics

Levey-Jennings charts along with cumulative mean and CVs are automatically calculated for each analyte and displayed on screen.

Supplemental Quality Monitoring (SQM)

Stat Profile Prime continuously monitors the status and performance of all analytical components (including sensors, reagents, calibrations, sample integrity, software, and electronics) providing real-time, sample-to-sample assurance of correct performance.

Compact Size

Stat Profile Prime has a 30.5 cm (12 in.) x 36.20 cm (14.35 in.) footprint to occupy minimal bench space. It is so compact it can also be located on a mobile cart for testing outside the lab.

Onboard Printer

Patient and QC results can be printed automatically or on demand. The printout also displays normal ranges for each test and flags results that exceed the normal range.

Plug-and-Play Connectivity to LIS/HIS

Stat Profile Prime easily connects to laboratory information systems and hospital information systems via industry standard, unidirectional POCT1-A2 or ASTM communication formats.



Choice of Stat Profile Prime® Models

Five models offer a choice of test menus from electrolytes only to a comprehensive, 10-test blood gas, electrolyte, and metabolite menu.

Prime Electrolyte System



Comprehensive Electrolyte Model*

Na, K, Cl, iCa, iMg

Basic Electrolyte Model*

Na, K, Cl

Acceptable Samples

Whole blood (heparinized), serum, plasma, urine, CSF (Chloride)

Sample Volume

Comprehensive Electrolyte Model 100 µL
Basic Electrolyte Model..... 50 µL

Calculated Parameters*

niCa, niMg, niCa/niMg

Auto-Sampler: An optional, factory-installed Auto-Sampler is available for the Stat Profile Prime electrolyte system. The ten-position tray accommodates serum, plasma, and urine samples in 2.0 ml and 0.5 ml sample cups. The Auto-Sampler is easy to load with prompting by the Stat Profile Prime user interface.

* Not yet available in the U.S.

Other Prime Critical Care Systems



Blood Gas/Electrolyte/Metabolite Model

pH, PCO₂, PO₂, Hct, Na, K, Cl, iCa, Glu, Lac

Blood Gas/Electrolyte Model

pH, PCO₂, PO₂, Hct, Na, K, Cl, iCa

Blood Gas Model

pH, PCO₂, PO₂

Acceptable Samples

Whole blood (heparinized), arterial, mixed venous, capillary

Sample Volume

Blood Gas Model 50 µL
Blood Gas/Electrolyte Model 100 µL
Blood Gas/Electrolyte/
Metabolite Model 100 µL

Calculated Parameters

SO₂%, HCO₃⁻, TCO₂, Be-efc, Be-b, SBC, O₂Ct, O₂Cap, A, AaDO₂, a/A, RI, PO₂/FIO₂, Anion Gap**, P50**, Hb**
Temperature corrected pH, PCO₂, PO₂

**Not available on the blood gas model

Measurement Range

Na	80 - 200 mmol/L	Glu	15 - 500 mg/dL (0.83 - 27.75 mmol/L)	pH	6.50 - 8.00
K	1 - 20 mmol/L	Lac	0.4 - 20 mmol/L (2.7 - 178.1 mg/dL)	PCO ₂	3.0 - 200 mmHg (0.4 - 26.6 kPa)
Cl	50 - 200 mmol/L	Hct	12% - 70%	PO ₂	5 - 765 mmHg (0.0 - 106.4 kPa)
iCa	0.2 - 2.7 mmol/L (0.8 - 10.3 mg/dL)			BarP	400.0 - 800.0 mmHg (53.3 - 106.7 kPa)
iMg	0.1 - 1.5 mmol/L				

Certifications ISO 9001 Quality System Registration, CSA, TÜV, CE Self Declared Complies to EN 61010, EN 50081, 82

Physical Specifications

Height: 39.06 cm (15.38 in.) Width: 30.5 cm (12 in.) Depth: 36.20 cm (14.35 in.) Weight: 17.9 lb (8.167 kg) without calibration cartridge

Optional Accessories

Combined 1D/2D Barcode Scanner: An optional, factory-installed barcode scanner is available for all Stat Profile Prime models. The scanner reads both 1D and 2D barcodes for patient and operator IDs. QC package inserts can also be scanned for lot number and expiration date.

nova
biomedical

NOVA BIOMEDICAL HEADQUARTERS: 200 Prospect Street, Waltham, MA 02453-3465 U.S.A. TEL: (781) 894-0800 (800) 458-5813 FAX: (781) 894-5915 Int'l FAX: (781) 899-0417 e-mail: info@novabio.com
NOVA BIOMEDICAL BRASIL: Avenida Barão Homem de Melo, 4391. Sala 608, 30494-275 Bairro Estoril, Belo Horizonte – MG – Brasil TEL: (55) 31 3657-0596 FAX: (55) 31 3342-2549 e-mail: novabrasil@novabiomedical.com.br
NOVA BIOMEDICAL CANADA, LTD: 17 – 2900 Argentia Road, Mississauga, Ontario L5N 7X9 Canada TEL: (905) 567-7700 (800) 263-5999 FAX: (905) 567-5496 e-mail: info@novabio.ca
NOVA BIOMEDICAL FRANCE: Parc Technopolis – Bât. Sigma 3 Avenue du Canada 91940 Les Ulis Courtaboeuf, France TEL: (33) 1-64 86 11 74 FAX: (33) 1-64 46 24 03 e-mail: info@novabiomedical.fr
NOVA BIOMEDICAL GmbH, Deutschland: Messenhäuser Str. 42, 63322 Rödermark, Germany TEL: (49) 6074-8448-0 FAX: (49) 6074-844833 e-mail: info@novabiomedical.de
NOVA BIOMEDICAL K.K.: Mita 43MT Building-7F, 13-16 Mita 3-chome, Minato-ku, Tokyo 108-0073, Japan TEL: (81) 3-5418-4141 FAX: (81) 3-5418-4676 e-mail: novabio@novabiomedical.co.jp
NOVA BIOMEDICAL U.K.: Innovation House, Aston Lane South, Runcorn, Cheshire WA7 3FY United Kingdom TEL: (44) 1928 704040 FAX: (44) 1928 796792 e-mail: sales@novabiomedical.co.uk
U.S.A. (Toll-Free) 800-458-5813 e-mail: info@novabio.com
CANADA (Toll-Free) 800-263-5999 www.novabiomedical.com